



ADVENTURER™

Semi-Micro, Analytical, and Precision Balances



By Model

Intuitive Balances Designed for Routine Weighing

Striking the ideal balance between inventive features and functional, uncomplicated weighing capabilities, the OHAUS Adventurer incorporates all of the applications necessary for routine weighing and measurement activities. With a color touchscreen, three level user management to fulfill GLP/GMP compliance capabilities, two USB ports, and much more, Adventurer is the most complete balance in its class.

Unique Features Include:

- Adventurer balances feature a color touchscreen, icon-based user interface, and an ergonomic design -making them easy to configure and use.
- Features such as specialized weighing modes, multiple connectivity options, and AutoCal™ provide versatility and flexibility for a variety of applications.
- Durable construction, large weighing surfaces, a space-saving draftshield design, and full housing in-use cover allow for use in lab, education and industrial environments.

Distributed by:
Safety Emporium
PO Box 1003
Blackwood, NJ 08012
Ph: (866) 326-5412 toll-free
Fax: (856) 553-6154
esupport@safetyemporium.com
www.safetyemporium.com

ADVENTURER™ *Semi-Micro, Analytical, and Precision Balances*

Stability, Accuracy, and Fast Operation Ensure Optimal Weighing Results in Routine Weighing Tasks

Weighing Performance

- Delivers stable and reliable weighing results for routine weighing tasks

Stabilization Time

- Adventurer's fast stabilization time improves productivity in the laboratory

Calibration

- AutoCal™ — Selected models feature OHAUS' automatic internal calibration system that performs routine maintenance by calibrating the balance daily
- External Calibration — Traditional calibration in which the operator manually calibrates the balance with their choice of calibration weight value to ensure accuracy available on every model



Color Touchscreen Offers Easy and Fast Operation of Adventurer's Applications

- Operate and access Adventurer's nine application modes and abundant features that eliminate the need to do several manual calculations through the modern color touchscreen
- Operators can wear laboratory gloves while utilizing the touchscreen, eliminating the inconvenience and hazards associated with constantly putting on and removing gloves
- In addition to the touchscreen, Adventurer also has six mechanical keys that provide tactile feedback and allow the operator to perform repetitive operations such as tare, zero, calibration, and print



Distributed by:
Safety Emporium
PO Box 1003
Blackwood, NJ 08012
Ph: (866) 326-5412 toll-free
Fax: (856) 553-6154
esupport@safetyemporium.com
www.safetyemporium.com

Application Modes



Weighing

Determine the weight of items in the selected unit of measure.



Parts Counting

Count samples of uniform weight.



Percent Weighing

Measure the weight of a sample displayed as a percentage of a pre-established Reference Weight.



Dynamic Weighing

Weigh an unstable load. Scale takes an average of weights over a period of time.



Density Determination

Determine density of solids or liquid. With the weigh below hook, it's possible to perform specific gravity tests for objects that cannot be easily placed on the weighing pan.



Check Weighing

Compare the weight of a sample against target limits.



Display Hold*

Manually holds the last stable weight or highest weighing value on the display.



Totalization / Statistics

Measure cumulative weight of multiple items. Cumulative total may exceed balance capacity.



Formulation

For compounding and recipe making. The number of components can range from 2 to 50.

*NTEP models (AXxxxN) will not have display hold.



Batch Printing

Combine multiple samples into one printout rather than printing them one at a time.

ADVENTURER™ *Semi-Micro, Analytical, and Precision Balances*

Equipped with the Connectivity and Functional Features Required in Laboratories

Distributed by:
Safety Emporium
PO Box 1003
Blackwood, NJ 08012
Ph: (866) 326-5412 toll-free
Fax: (856) 553-6154
esupport@safetyemporium.com
www.safetyemporium.com

Dual USB Ports

- A front USB host port is easily accessible and makes it simple to load data from the balance on to a flash drive without having to reach around to the back or move the balance
- A second USB device port is located at the rear of the balance that can be used to connect the balance to a PC
- The connectivity options help meet traceability requirements in traditional installations

Label Printing Function

- Easy to link with Zebra printer and available with one built-in label printing template

Balance Profiles

- The cloning feature allows you to save user and application settings to a USB flash drive which can be easily used to configure additional Adventurer balances

Below Minimum Sample Weight Indication

- When using the minimum weight feature, the display clearly indicates that your current sample weight is below your defined minimum limit. Simply increase your sample weight to assure that your results are up to your standard



Space-saving Draftshield Designed to Improve User Experience and Accessibility

- Draftshield doors are constructed of two glass panels, reducing the space required on the lab bench when the doors are open
- Wide door entry provides unobstructed access and allows larger weighing vessels to be easily placed on the pan, reducing the chance of accidental spillage
- Easy to keep clean in order to ensure a safe workspace by minimizing contamination



Power Saving Functions

- The Adventurer is designed with various power-saving features that help reduce the environmental impact. These include auto-off and brightness controls.



User Management Function ensures data security and data traceability

- 3 level user management function ensures data security and data traceability requirements
- One administrator, two supervisors and 10 users have preset accessibility in the software



Real Time Clock with GLP/GMP Data

- A real-time clock function keeps accurate time even during power loss
- GLP data capability has the ability to record Sample name, Project names and Balance IDs to help meet traceability and compliance requirements

ADVENTURER™ *Semi-Micro, Analytical, and Precision Balances*

Specifications

Model	AX125D	AX85	AX225D	AX124	AX224	AX324	AX223	AX423	AX523	AX623
External Calibration	—	—	—	AX124/E	AX224/E	—	AX223/E	AX423/E	AX523/E	AX623/E
Approved Models	—	—	—	—	AX224N	—	AX223N/E	AX423N AX423N/E	AX523N/E	AX623N/E
Capacity (g)	82/120	82	102/220	120	220	320	220	420	520	620
Readability d (g)	0.00001/0.0001	0.00001	0.00001/0.0001	0.0001			0.001			
Verification Interval* e (g)	—			—	0.001	—	0.01			
Class*	—			I			II			
*Repeatability (sd.), ≤5% of Full Load (g)"	0.00001			0.00008			0.0008			
Repeatability (sd.), 5% of Full Load to Fine Range Max (g)	0.00002			—						
Repeatability (sd.), 5% of Full Load to Full Range (g)	0.0001	0.00002	0.0001	0.0001			0.001			
Linearity Deviation, Typical (g)	±0.00006	±0.00006	±0.00006	±0.00006			±0.0006			
Linearity Deviation (g)	±0.0001	±0.0001	±0.0001	±0.0002			±0.002			
Stabilization Time (sec)	8s			≤3			≤2			
Sensitivity Drift (ppm/°C)	2			3			3			
Min-Weight (Typical) (g) (USP, K=2, U=0.10%)	20 mg			0.16			0.82 g			
*Min-Weight (Optimal) (g) (USP, K=2, U=0.10%, SRP≤0.41d)***	8.2 mg			0.082 g			1.6			
Weighing Units	gram, milligram, kilogram, mesgal, momme, Newton, ounce, pennyweight, Baht, carat, grain, pound, Tael (Hong Kong), Tael (Singapore), Tael (Taiwan), tical, tola, troy ounce, custom (1)									
Weighing Units, Approved Models	—			ct, grain, g, mg, oz, ozt, dwt						
Weighing Applications	Weighing, Parts Counting, Percent Weighing, Check Weighing, Dynamic Weighing, Formulation, Density Determination, Totalization, Display Hold, Batch Printing									
Pan Size	Ø 3.1 in / 80 mm			Ø 3.5 in / 90 mm			Ø 5.1 in / 130 mm			
Calibration	All models feature external calibration. Models feature AutoCal™ internal calibration, except for AX.../E models									
Tare Range	To capacity by subtraction									
Power Requirements	AC Adapter Input: 100-240 VAC 0.3A 50-60 Hz AC Adapter Output: 12 VDC 0.84A									
Display Type	Full-Color Touchscreen WQVGA Graphic LCD									
Display Size	4.3 in / 109 mm (diagonal)									
Base Housing (W×H×D)	13.93 x 13.39 x 9.05 in / 354 x 340 x 230 mm									
Communication	RS232, USB Device, USB Host									
Temperature Range	50°F to 86°F/10°C to 30°C									
Humidity Range	Maximum relative humidity 80 % for temperatures up to 86°F/30°C									
Storage Conditions	14°F/-10°C to 140°F/60°C at 10% to 90% relative humidity, non-condensing									
Shipping Dimensions	20.0 x 15.4 x 20.9 in / 507 x 387 x 531 mm									
Net Weight	11.3 lb / 5.1 kg						12.8 lb / 5.8 kg			
Shipping Weight	17.2 lb / 7.8 kg						18.8 lb / 8.5 kg			

Model	—	AX622	—	AX1502	AX2202	AX4202	AX5202	AX6202	—	AX4201	AX8201	AX12001	
External Calibration	AX422/E	—	AX822/E	AX1502/E	AX2202/E	AX4202/E	—	AX6202/E	AX2201/E	AX4201/E	AX8201/E	AX12001/E	
Approved Models	—	AX622N/E	—	AX1502N/E	AX2202N/E	AX4202N/E	—	AX6202N/E	—	—	AX8201N/E	AX12001N/E	
Capacity (g)	420	620	820	1,520	2,200	4,200	5,200	6,200	2,200	4,200	8,200	12,000	
Readability d (g)	0.01												
Verification Interval* e (g)	—				0.1				—				1
Class*	II			—				II		—			II
*Repeatability (sd.), ≤5% of Full Load (g)"	0.008												
Repeatability (sd.), 5% of Full Load to Fine Range Max (g)	—												
Repeatability (sd.), 5% of Full Load to Full Range (g)	0.01						0.1						
Linearity Deviation, Typical (g)	±0.006						±0.06						
Linearity Deviation (g)	±0.02						±0.2						
Stabilization Time (sec)	≤1.5												
Sensitivity Drift (ppm/°C)	3						5						
Min-Weight (Typical) (g) (USP, K=2, U=0.10%)	16						160						
*Min-Weight (Optimal) (g) (USP, K=2, U=0.10%, SRP≤0.41d)***	8.2 g						82 g						
Weighing Units	gram, kilogram, mesgal, momme, Newton, ounce, pennyweight, Baht, carat, grain, pound, Tael (Hong Kong), Tael (Singapore), Tael (Taiwan), tical, tola, troy ounce, custom (1)												
Weighing Units, Approved Models	ct, grain, g, oz, ozt, dwt, lb, kg												
Weighing Applications	Weighing, Parts Counting, Percent Weighing, Check Weighing, Dynamic Weighing, Formulation, Density Determination, Totalization, Display Hold, Batch Printing												
Pan Size	6.9 x 7.7 in / 175 x 195 mm												
Calibration	All models feature external calibration. Models feature AutoCal™ internal calibration, except for AX.../E models												
Tare Range	To capacity by subtraction												
Power Requirements	AC Adapter Input: 100-240 VAC 0.3A 50-60 Hz AC Adapter Output: 12 VDC 0.84A												
Display Type	Full-Color Touchscreen WQVGA Graphic LCD												
Display Size	4.3 in / 109 mm (diagonal)												
Base Housing (W×H×D)	13.93 x 3.94 x 9.05 in / 354 x 100 x 230 mm												
Communication	RS232, USB Device, USB Host												
Temperature Range	50°F to 86°F/10°C to 30°C												
Humidity Range	Maximum relative humidity 80 % for temperatures up to 86°F/30°C												
Storage Conditions	14°F/-10°C to 140°F/60°C at 10% to 90% relative humidity, non-condensing												
Shipping Dimensions	22.0 x 15.5 x 11.9 in / 557 x 392 x 301 mm												
Net Weight	11.3 lb / 5.1 kg						12.8 lb / 5.8 kg						
Shipping Weight	17.2 lb / 7.8 kg						18.8 lb / 8.5 kg						


Distributed by:
Safety Emporium
 PO Box 1003
 Blackwood, NJ 08012
 Ph: (866) 326-5412 toll-free
 Fax: (856) 553-6154
esupport@safetyemporium.com
www.safetyemporium.com

*Certified models only **The value for SRP is the standard deviation for n replicate weighings (n≥10)


ADVENTURER™ *Semi-Micro, Analytical, and Precision Balances*

Outline Dimensions

0.1 and 1mg Models

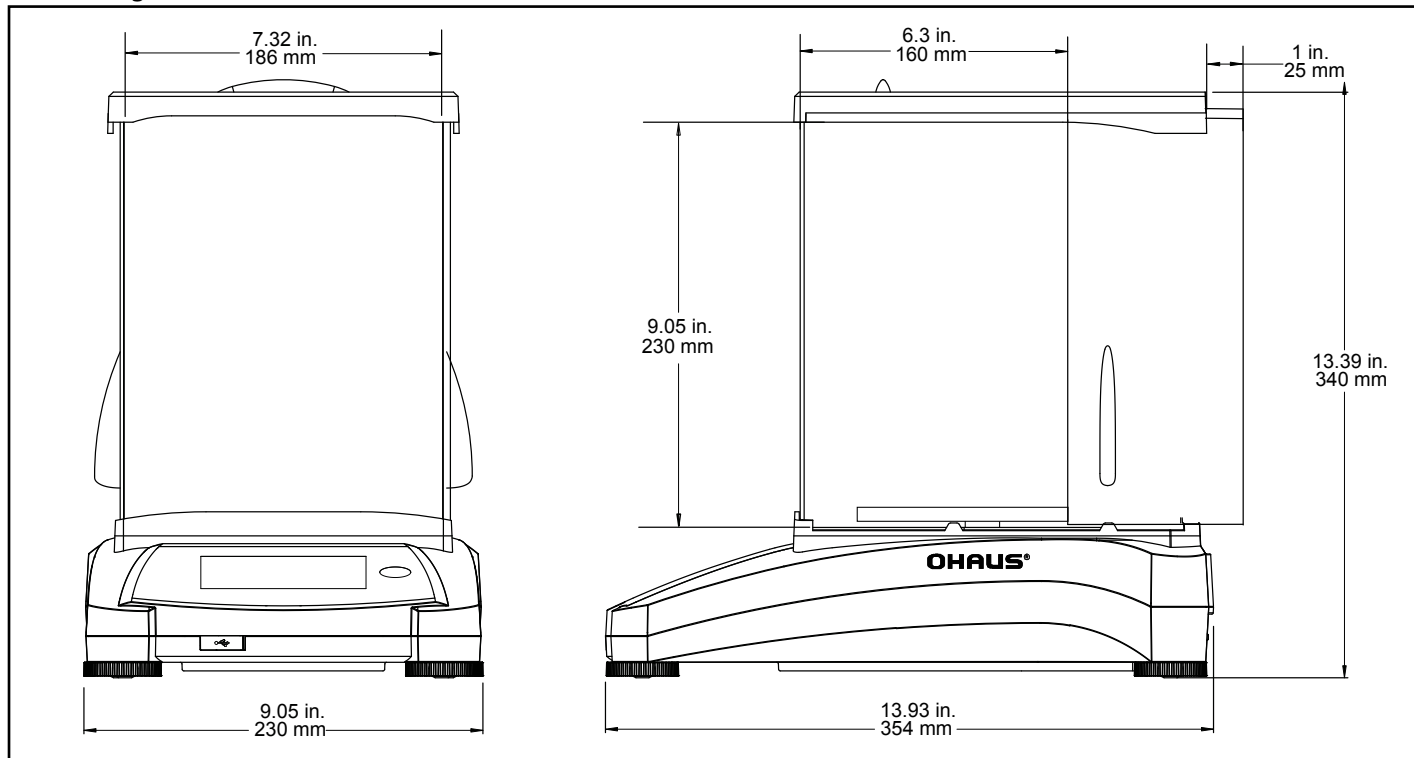

0.01 and 0.1g Models

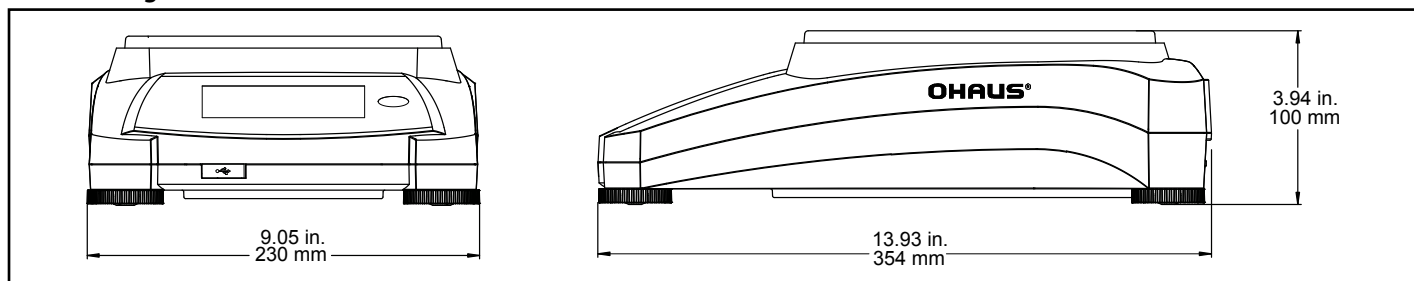

Additional Features

RS232 interface, integrated weigh below hook, full housing in-use cover, removable stainless steel pan, die-cast metal bottom housing, security bracket, illuminated up-front level indicator, four adjustable feet, software lockout menus, stability indicator, software overload/underload indicators, user selectable environmental settings, audible indicator, user selectable brightness settings, auto dim, auto-standby, auto-off, touchscreen calibration, auto tare, user selectable operating language (14), compatible interface command with MT-SICS and ST protocol

Compliance

Metrology: NIST Handbook 44, Measurement Canada Weights and Measures Regulations (Class I, nmax 220000; Class II, nmax 62000)

Product Safety: CAN/CSA C22.2 61010-1, UL 61010-1, IEC 61010-1

Electromagnetic Compatibility: FCC Part 15 Class A, ICES-001 Class A, IEC 61326-1 (emissions Class B, immunity Basic requirements)

Accessories

ION-100A US Standalone Ionizer	30130302
STP103 Printer	80251992
SF40A Impact Printer	30064203
Auxiliary Display	30472064
Density Determination Kit	80253384
Cable, USB Interface (Type A to B)	83021085
Security Device (Laptop Lock)	80850043
RS232 Cable, PC 9 Pin	80500525

Distributed by:
 Safety Emporium
 PO Box 1003
 Blackwood, NJ 08012
 Ph: (866) 326-5412 toll-free
 Fax: (856) 553-6154
 esupport@safetyemporium.com
 www.safetyemporium.com

80774742_M 20241113 © Copyright OHAUS Corporation

OHAUS CORPORATION

Tel: 800.672.7722
 973.377.9000
 Fax: 973.944.7177

www.ohaus.com

The management system governing the manufacture of this product is ISO 9001:2015 certified.

